



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Winston-Lutz Analysis

Number of images: 9

Maximum 2D CAX->BB distance: 1.42mm

Median 2D CAX->BB distance: 0.93mm

Shift to iso: facing gantry, move BB: RIGHT 0.33mm; IN 0.48mm; DOWN 0.20mm

Gantry 3D isocenter diameter: 1.18mm (4/9 images considered)

Maximum Gantry RMS deviation (mm): 0.73mm

Maximum EPID RMS deviation (mm): 2.30mm








Gantry+Collimator 3D isocenter diameter: 1.60mm (7/9 images considered)

Collimator 2D isocenter diameter: 0.61mm (4/9 images considered)

Maximum Collimator RMS deviation (mm): 1.42

Couch 2D isocenter diameter: 1.28mm (3/9 images considered)

Maximum Couch RMS deviation (mm): 1.34

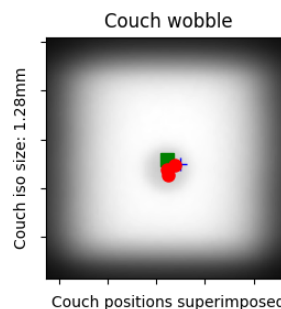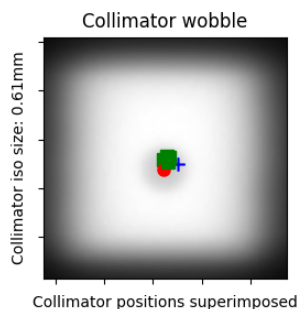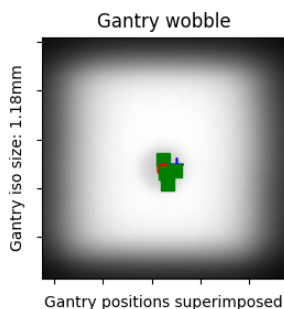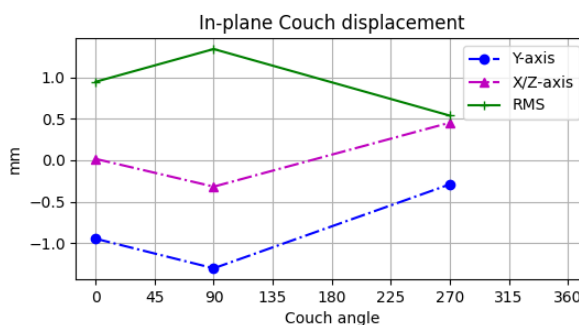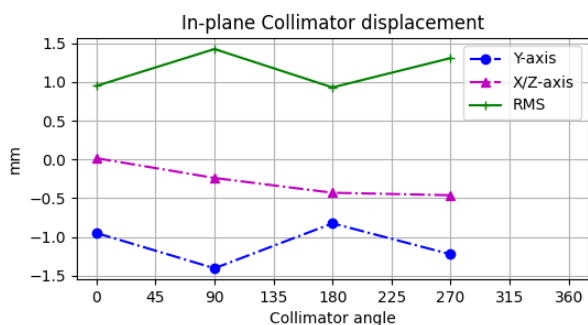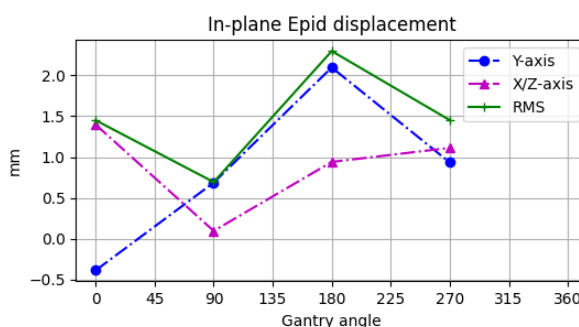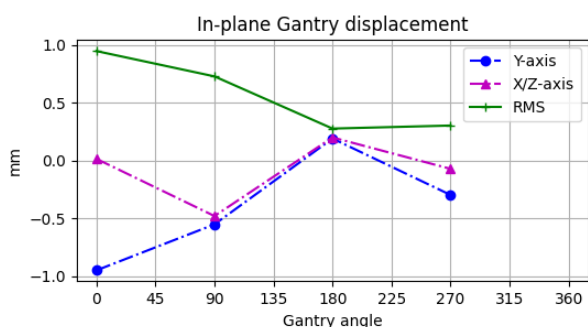

Notes:

Análise de Winston-Lutz com pylinac

Campo 2 cm X 2 cm

Tolerância: 1 mm

Limite de ação: 0.8 mm

Local das imagens: ./WL_IOSP\WL_2021-05-11a



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Gantry images

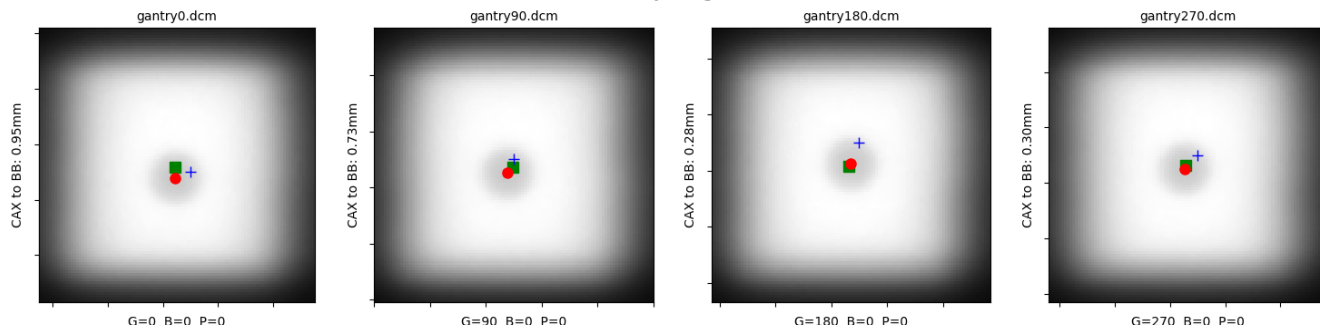



Winston-Lutz Analysis

Metadata:


Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Collimator images

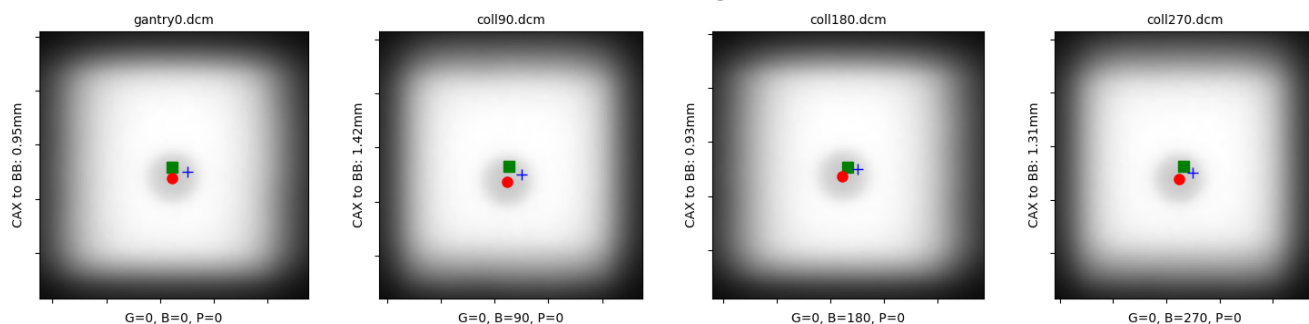



Winston-Lutz Analysis

Metadata:

Autor: William A. P. dos Santos

Função: Físico

Acelerador Linear: Synergy Full

Modelo MLC: Agility

Couch images

